





PET TOOTH CLEANER

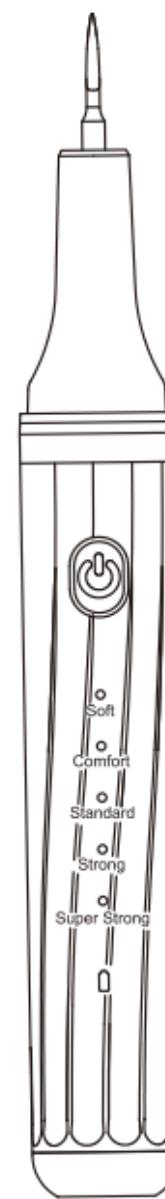

1. Notice

- (1) This product uses ultrasonic technology. There is no vibration or sound when turned on; it only works upon contact with teeth or hard objects and stops automatically when removed.(If it doesn't work, please turn off the device, tighten the cleaning head with the wrench, and turn it back on.)
- (2) The product has an automatic shutdown feature. It turns off after 9 minutes to prevent overheating, extending product life. You can restart it at any time after shutdown.
- (3) Adjust usage time according to your pet's dental condition, as frequent use may damage their teeth.
- (4) During use, the product may heat up (cleaning head up to 52°C, handle up to 43°C), which is normal. The cleaning head cools down in 10 seconds after stopping, and the handle cools down in 8 minutes.

2. Packing List

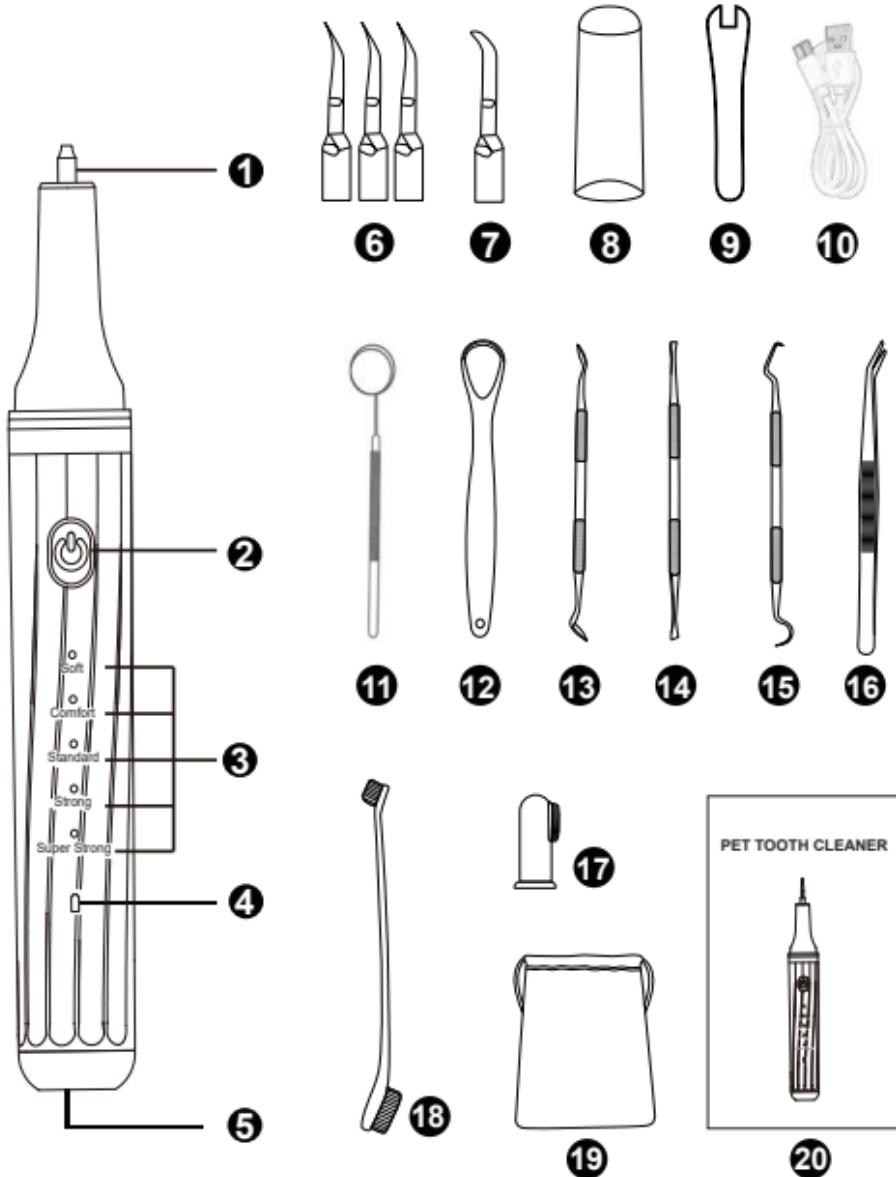

- | | |
|---------------------------|----------------------------------|
| ① LED Light | ⑪ 1*Oral Mirror |
| ② Power Button | ⑫ 1*Tongue Scraper |
| ③ Mode Indicator | ⑬ 1*Sickle Type Tooth
Cleaner |
| ④ Power Indicator | ⑭ 1*Hoe Type Tooth Cleaner |
| ⑤ USB Charging Port | ⑮ 1*Probe Teeth Cleaner |
| ⑥ 3*Pointed Cleaning Head | ⑯ 1*Tweezers |
| ⑦ 1*Flat Cleaning Head | ⑰ 1*Pet Finger Toothbrush |
| ⑧ Dust Cover | ⑱ 1*Pet Two-headed
Toothbrush |
| ⑨ 1*Wrench | ⑲ 1*Flannel Storage Bag |
| ⑩ 1*USB Type-C Charging | ⑳ 1*Specification |

3. Charge Before Use

Due to long-distance transportation, the power will be reduced. Please fully charge the device before use
(Fully charged in 2 hours)

- (1) The power indicator light is red when charging.
- (2) The power indicator light is green after fully charged.
- (3) During the charging, the device cannot work.
- (4) Fully charge in 2 hours. It is recommended to use a 5V 1A adapters.

4. Operation Guide

- (1) Power on the Device: Press and hold the power button to turn on the device. The indicator light will turn white, and the LED lights will turn on. Press the power button to switch modes. Choose from five modes: soft, comfort, standard, strong, and super strong, depending on your pet's dental condition.
- (2) Power off: After use, press and hold the power button for 2-3 seconds; the indicator light will go out, and the device will turn off.
- (3) During use, if the power indicator turns red, the battery is low. Please charge it soon. If the red light flashes, the device will stop working and needs charging.

5.Cleaning Head Introduction


- (1) Pointed Cleaning Head: Used to remove calculus and tooth smoke stains on gingiva.
- (2) Flat Cleaning Head: Used to remove large stones and stains on the teeth and lingual side of the supra gingival teeth.

6.How to Use Cleaning Head


Gently press the cleaning head vertically against the teeth and move it up and down, as well as left and right. Adjust the usage time for each area based on the condition of the teeth.

Note: Please do not use this product on teeth without tartar.


(1) Tooth surface dirt

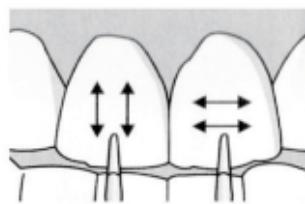


(2) Boundary between gum and tooth

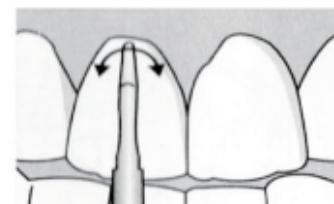


(3) Interdental gap

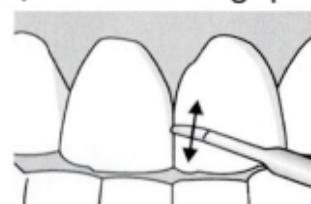


(4) Dirt on uneven teeth

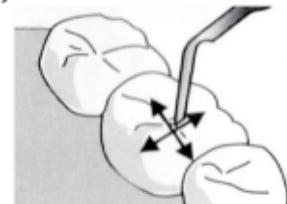

(5) Dirt on the back of teeth

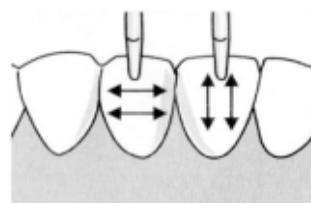

(6) Use the oral mirror when cleaning the back of teeth

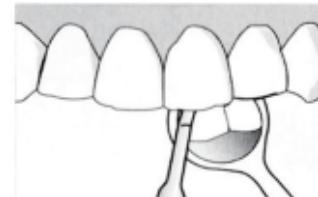

(7) Use the tongue scraper to clean your pet's tongue coating



(8) Use tweezers to pick up waste after cleaning




- (9) Use a probe to clean between your pet's teeth (10) Use a hoe-type tooth cleaner for surface stones and stains


- (11) Use a sickle-type tooth cleaner for stones under the gums

7. How to Install the Cleaning Head

The orientation of the cleaning head may be different, but as long as it is firmly mounted on it, it can be used in any direction.

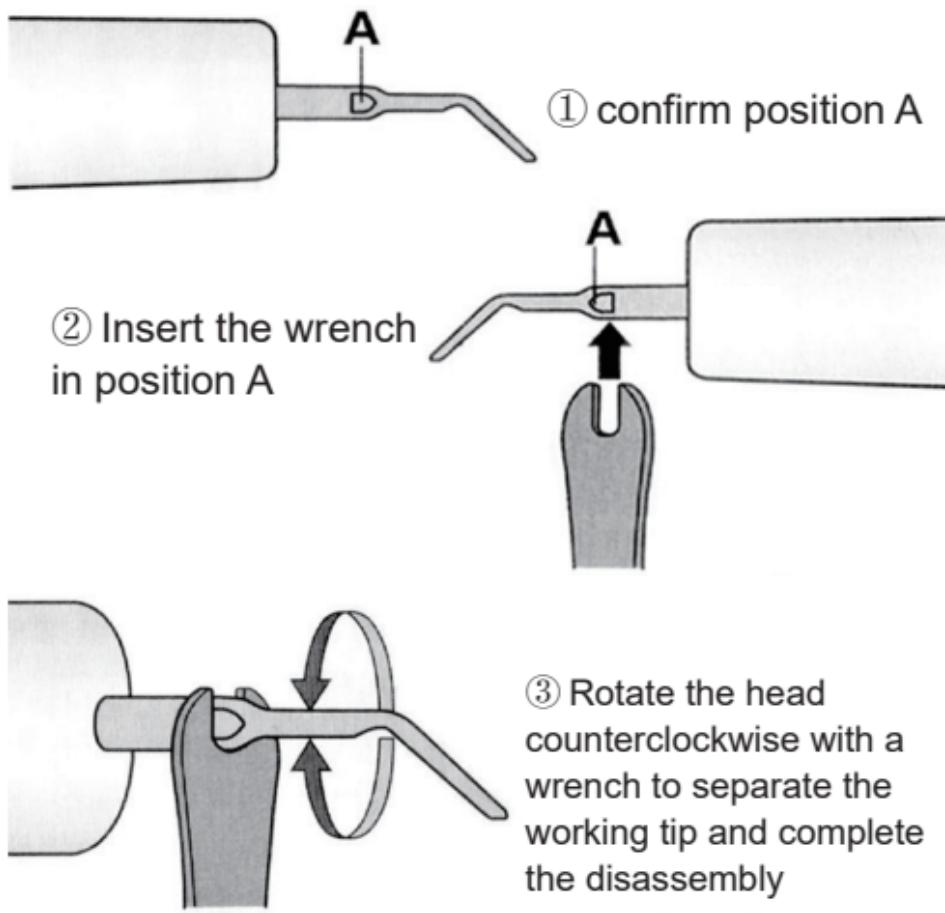

Note: Preventive measure
Avoid over tight assembly resulting in damage to the interface.

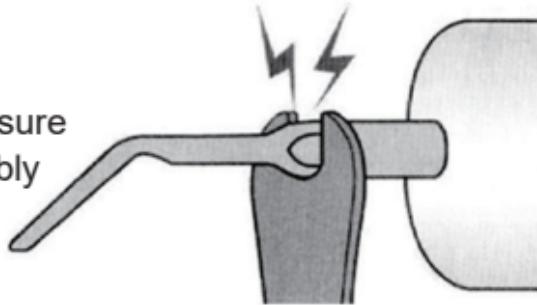

8.How to Clean

- (1) Before the cleaning head is used,it can be disinfected by steaming, alcohol soaking,high temperature and high pressure, etc.
- (2) When not in use for a long time, please keep it in a dry place.

	
The device can be wiped with rag	Can't be soaked in water for a long time

9.Specification

Charging Voltage	5V
Charging Current	≤1A
Charging Time	2Hours
Frequency	32KHZ-42KHZ
Consumption	≤2.5w
Battery Capacity	500mAh

Product name: Pet Tooth Cleaner

SKU: AEPETC0101/AEPETC0102

Input voltage/current: 5V/1A

Output power: 2.5W

Model: HS-03

⚠ Warning: Please use with care, keep away from pets and children, clean and maintain regularly. This product is for pets only. If any abnormal conditions occur, please stop using and contact a veterinarian.

Made in China

Manufacturer: Shenzhen Oulepu Technology Co., Ltd.

Manufacturing location: Shenzhen, China

Address: Room 301, Yuanxiong Building, Bantian Dadao, Bantian street, Longgang District, Shenzhen

Manufacturer E-mail:

customer-support@petbank-co.com

[EC REP]

EU representative: AMZLAB GmbH

Address: Laubenhof 23, 45326 Essen

Customer Service Email:

customer-support@petbank-co.com

[UK AR]

UK representative: UK CROSSBORDER LIMITED

Address: 7 Bell Yard London WC2A 2JR, UK United Kingdom

Customer Service Email:

customer-support@petbank-co.com

